

PERIPHERALLY INSERTED CENTRAL CATHETER (PICC)

## VASCU-PICC® 2.6F

### 1.9F SINGLE LUMEN PICC

Catalog Code	Description	QTY
MR17011100	1.9F x 50CM SINGLE LUMEN PICC, CATH ONLY	10/BOX

#### Tray Includes

- |                                   |                                 |
|-----------------------------------|---------------------------------|
| (1) Vascu-Picc W/Preloaded Stylet | (1) Patient Information Booklet |
| (1) Needle Free Valve             | (1) Patient Chart Sticker       |
| (1) Catheter Securement Device    |                                 |

Catalog Code	Description	QTY
MR17011121	1.9F x 20CM SINGLE LUMEN PICC	5/BOX
MR17011101	1.9F x 50CM SINGLE LUMEN PICC	5/BOX

#### Tray Includes

- |                                   |                                |
|-----------------------------------|--------------------------------|
| (1) Vascu-Picc W/Preloaded Stylet | (1) 2F (22ga) Tearaway         |
| (1) Latex Free Tourniquet         | (1) 24" Tape Measure           |
| (10) 2"x2" Gauze                  | (1) Catheter Securement Device |
| (1) 10CC Syringe                  | (1) Patient Chart Sticker      |
| (1) Needle Free Valve             |                                |

### 2.6F DUAL LUMEN PICC

Catalog Code	Description	QTY
MR17012600	2.6F x 50CM DOUBLE LUMEN PICC, CATH ONLY	10/BOX

#### Tray Includes

- |                                   |                                 |
|-----------------------------------|---------------------------------|
| (1) Vascu-Picc W/Preloaded Stylet | (1) Patient Information Booklet |
| (2) Needle Free Valve             | (1) Patient Chart Sticker       |
| (1) Catheter Securement Device    |                                 |

Catalog Code	Description	QTY
MR17012621	2.6F x 20CM DOUBLE LUMEN PICC	5/BOX
MR17012601	2.6F x 50CM DOUBLE LUMEN PICC	5/BOX

#### Tray Includes

- |                                   |                                |
|-----------------------------------|--------------------------------|
| (1) Vascu-Picc W/Preloaded Stylet | (1) 3F (19ga) Tearaway         |
| (1) Latex Free Tourniquet         | (1) 24" Tape Measure           |
| (10) 2"x2" Gauze                  | (1) Catheter Securement Device |
| (1) 10CC Syringe                  | (1) Patient Chart Sticker      |
| (2) Needle Free Valves            |                                |

### PROCEDURE TRAY

Catalog Code	Description	QTY
MR17019000	PROCEDURE TRAY	3/BOX

#### Tray Includes

- |                                      |                                 |                               |
|--------------------------------------|---------------------------------|-------------------------------|
| (3) Pre-Filled Saline Syringes       | (1) #11 Safety Scalpel          | (1) Scissors                  |
| (2) Chlorhexidine                    | (1) Plastic Forceps             | (1) Surgical Tape             |
| (1) 25" x 14" Absorbant Towel        | (1) 37" x 24" Fenestrated Drape | (1) Steri-Strips              |
| (1) 36" x 36" CSR Wrap               | (1) Adhesive Dressing           | (1) Needle Stick Pad          |
| (1) 5CC Lidocaine HCL 1%             | (1) Filter Straw                | (1) Lidocaine/Saline Stickers |
| (3) Surgical Gloves                  | (1) Saline Wipe                 | (1) 24" Tape Measure          |
| (1) 40" x 70" Drape                  | (2) Face Masks                  | (1) 36" x 24" Drape           |
| (1) 5CC Disposable Luer Lock Syringe | (1) Bouffant Cap                | (10) 4" x 4" Gauze            |
| (1) Hypodermic Safety Needle         | (1) 24" x 24" CSR Wrap          | (9) 2" x 2" Gauze             |
| (1) Gown                             |                                 |                               |

### PEDIATRIC PROCEDURE TRAY


Catalog Code	Description	QTY
MR17019001	PEDIATRIC PROCEDURE TRAY	3/BOX

#### Tray Includes

- |                                      |                                 |                               |
|--------------------------------------|---------------------------------|-------------------------------|
| (3) Pre-Filled Saline Syringes       | (1) Gown                        | (1) Scissors                  |
| (2) Chlorhexidine                    | (1) #11 Safety Scalpel          | (1) Surgical Tape             |
| (1) 25" x 14" Absorbant Towel        | (1) Steel Forceps               | (1) Steri-Strips              |
| (1) 36" x 36" CSR Wrap               | (1) 19" x 12" Fenestrated Drape | (1) Needle Stick Pad          |
| (1) 5CC Lidocaine HCL 1%             | (1) Adhesive Dressing           | (1) Lidocaine/Saline Stickers |
| (3) Surgical Gloves                  | (1) Filter Straw                | (1) 24" Tape Measure          |
| (1) 5CC Disposable Luer Lock Syringe | (1) Saline Wipe                 | (1) 36" x 24" Drape           |
| (1) Hypodermic Safety Needle         | (2) Face Masks                  | (10) 4" x 4" Gauze            |
| (1) 24" x 24" CSR Wrap               | (1) Bouffant Cap                | (9) 2" x 2" Gauze             |

 medCOMP®

*One size does not fit all.*

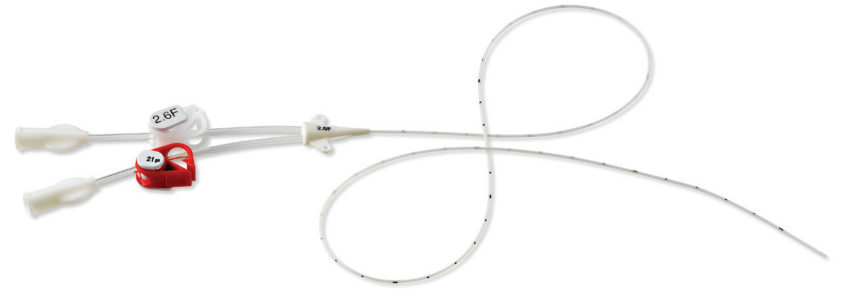

PERIPHERALLY INSERTED CENTRAL CATHETER (PICC)

## VASCU-PICC® 2.6F

## DUAL LUMEN VASCU-PICC® 2.6F

**SMALL but POWERFUL**

Approved for infusion therapy

**From the smallest to the largest — and now everyone in-between — Medcomp® has a catheter size for any pediatric patient.** VASCU-PICC® 2.6F is the latest addition to our PICC product line and features a dual lumen configuration, soft polyurethane body, and high flow rates in a small French size. Since smaller veins require greater precision, it only makes sense to offer you more choices in PICC sizes.

*Especially with your smallest patients.*

### Comfort

Comfort is a special challenge with your pediatric patients. The smaller French size features a dual lumen which means fewer sticks, fewer tears.

### Versatility

Since many pediatric patients require a combination of nutrition, fluids and/or medicine, the dual lumen configuration enables one catheter to serve multiple functions.

### Perfect Fit


Medcomp® offers a full-line of pediatric PICC sizes to fit any pediatric patient. The new 2.6F fits nicely between existing 1.9F and 3F catheters, and makes catheter insertion easier.

### 2.6F Catheter Size 20cm & 50cm Lengths

Smaller catheter while maintaining dual access and reduced trauma in smaller vessels

### Thermosensitive Material

Softens at body temperature

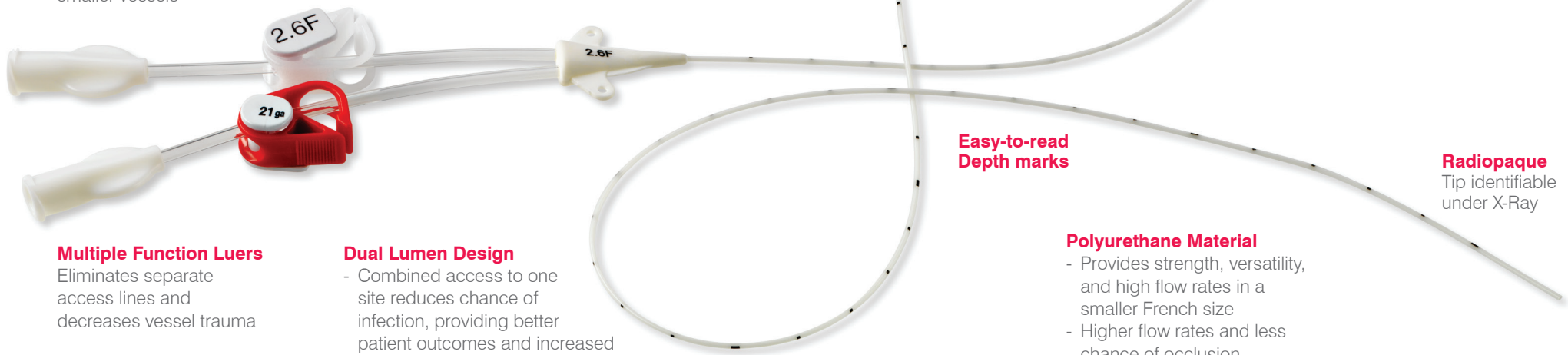

### Multiple Function Luers

Eliminates separate access lines and decreases vessel trauma

### Dual Lumen Design

- Combined access to one site reduces chance of infection, providing better patient outcomes and increased parent satisfaction
- Offers the ability to infuse incompatible drugs

### Easy-to-read Depth marks

### Polyurethane Material

- Provides strength, versatility, and high flow rates in a smaller French size
- Higher flow rates and less chance of occlusion
- Allows greater control during insertion

### Radiopaque Tip identifiable under X-Ray