




PRO-PICC<sup>®</sup>  
& VASCU-PICC<sup>®</sup>  
VALVED CATHETERS





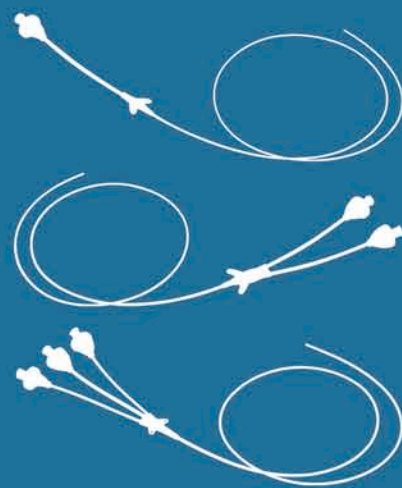


**Prevent  
catheter  
occlusion &  
reduce  
costs of care.**

The impact of catheter occlusions and infections can be measured not only by the quantifiable, but emotional costs that result from patient readmissions, additional days of treatment, catheter replacements, and patient discomfort and inconvenience. The incorporation of a proximal valve reduces the reflux of blood into the distal catheter tip after the device has been accessed or flushed, which can reduce the incidence of occlusion and infection. Although reductions in infection as great as 50% have been observed, a valved catheter becomes cost-effective with a reduction in occlusions of just 10% or infection of 1%.<sup>(1)</sup>

In addition, the costs of maintenance are reduced with the ability to eliminate heparin, using a saline-only lock instead. Medcomp's® unique valve design provides an effective solution, which optimizes the performance of the valve by reducing the interior low-flow region, which can be problematic with other designs.

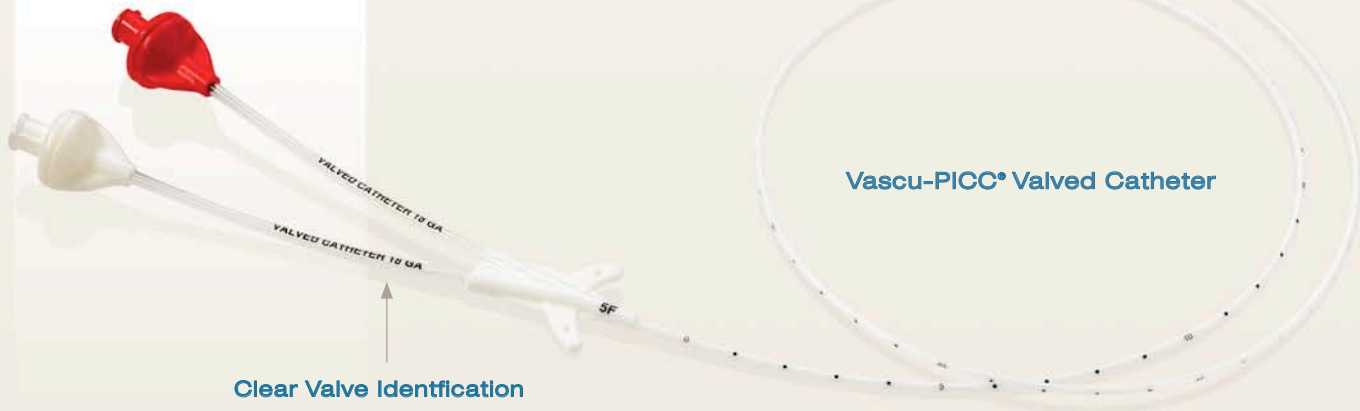

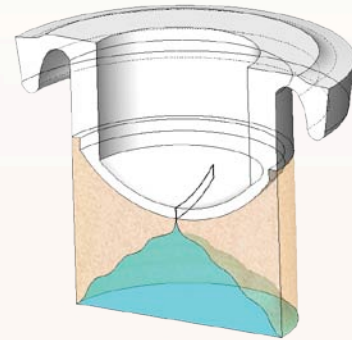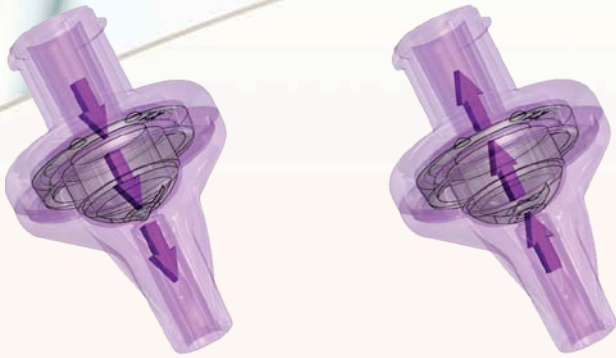

**Bi-directional performance** enables smooth passage of a wire in either direction as well as ease of aspiration and infusion.

**Minimal Low-Flow Regions.** Requires less negative pressure for easier blood return while also allowing for a more effective flush.

## FEATURES & BENEFITS

- **Clear Valve Identification**  
All valved PICCs are clearly labeled on the extensions.
- **CVP Monitoring**  
All valved PICCs are indicated for Central Venous Pressure Monitoring.
- **Multiple Kit Configurations**  
Available in a variety of kit configurations to accommodate insertion needs.
- **Saline Only Flush**  
Limits the potential complications & costs of heparin use.
- **Power Injection Rates**  
Achieve maximum injection rates of 5cc/sec at 300psi on all Pro-PICC Valved catheters.
- **Clamp Free Extensions**  
Eliminates maintenance challenges associated with needless injection caps.

## Pro-PICC® Valved Catheter

**Clamp Free Extensions**



**Power Injection Rates**

<sup>1</sup> Hoffer, MD, Eric K., John Borsa, et al. "Prospective Randomized Comparison of Valved Versus Nonvalved Peripherally Inserted Central Vein Catheters." AJR. 173. (1999): 1393-1398. Print.

## Pro-PICC® Valved Catheter

CATALOG #	CATHETER ONLY	UOM
MR82034100	4F X 55CM PRO-PICC® VALVED, SINGLE LUMEN	10/BOX
MR82035200	5F X 55CM PRO-PICC® VALVED, DOUBLE LUMEN	10/BOX
MR82036300	6F X 60CM PRO-PICC® VALVED, TRIPLE LUMEN	10/BOX

CATH ONLY CONTAINS: (1) PICC (1|2|3) Needle Free Valve(s) (1) Patient Information Pack (1) Patient Chart Sticker (1) Catheter Securement Device (1) Twisted Wire Stylet

CATALOG #	BASIC SET	UOM
MR82034101	4F X 55CM PRO-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82035101	5F X 60CM PRO-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82035201	5F X 55CM PRO-PICC® VALVED, DOUBLE LUMEN	5/BOX
MR82036301	6F X 60CM PRO-PICC® VALVED, TRIPLE LUMEN	5/BOX

BASIC SET CONTAINS: (1) PICC (1) Peelable Sheath Introducer (1) Mini Scalpel (1) Guidewire (1) 10cc Syringe (1) Tape Measure (1) Introducer Needle (1|2|3) Needle Free Valve(s) (1) Catheter Securement Device (1) Patient Information Pack (1) Patient Chart Sticker

CATALOG #	LONG WIRE SET	UOM
MR82034102	4F X 55CM PRO-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82035102	5F X 60CM PRO-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82035202	5F X 55CM PRO-PICC® VALVED, DOUBLE LUMEN	5/BOX
MR82036302	6F X 60CM PRO-PICC® VALVED, TRIPLE LUMEN	5/BOX

LONG WIRE SET CONTAINS: (1) PICC (1) Peelable Sheath Introducer (1) Mini Scalpel (1) Guidewire (1) 10cc Syringe (1) Tape Measure (1) Introducer Needle (1|2|3) Needle Free Valve(s) (1) Catheter Securement Device (1) Patient Information Pack (1) Patient Chart Sticker

CATALOG #	NURSING KIT	UOM
MR82034105	4F X 55CM PRO-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82035105	5F X 60CM PRO-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82035205	5F X 55CM PRO-PICC® VALVED, DOUBLE LUMEN	5/BOX
MR82036305	6F X 60CM PRO-PICC® VALVED, TRIPLE LUMEN	5/BOX

NURSING KIT CONTAINS: (1) PICC (1) Peelable Sheath Introducer (1) Safety Scalpel (1) Guidewire (1) 10cc Syringe (1) Tape Measure (1) Safety Introducer Needle (1|2|3) Needle Free Valve(s) (1) Catheter Securement Device (1) Latex-Free Tourniquet (1) Patient Information Pack (1) Patient Chart Sticker (1) Twisted Wire Stylet

## Vascu-PICC® Valved Catheter

CATALOG #	CATHETER ONLY	UOM
MR82013100	3F X 60CM VASCU-PICC® VALVED, SINGLE LUMEN	10/BOX
MR82014100	4F X 60CM VASCU-PICC® VALVED, SINGLE LUMEN	10/BOX
MR82014200	4F X 60CM VASCU-PICC® VALVED, DOUBLE LUMEN	10/BOX
MR82015200	5F X 60CM VASCU-PICC® VALVED, DOUBLE LUMEN	10/BOX

CATH ONLY CONTAINS: (1) PICC (1|2) Needle Free Valve(s) (1) Patient Information Pack (1) Patient Chart Sticker (1) Catheter Securement Device (1) Twisted Wire Stylet

CATALOG #	BASIC SET	UOM
MR82013101	3F X 60CM VASCU-PICC VALVED, SINGLE LUMEN	5/BOX
MR82014101	4F X 60CM VASCU-PICC VALVED, SINGLE LUMEN	5/BOX
MR82014201	4F X 60CM VASCU-PICC VALVED, DOUBLE LUMEN	5/BOX
MR82015201	5F X 60CM VASCU-PICC VALVED, DOUBLE LUMEN	5/BOX

BASIC SET CONTAINS: (1) PICC (1) Peelable Sheath Introducer (1) Mini Scalpel (1) Guidewire (1) 10cc Syringe (1) Tape Measure (1) Introducer Needle (1|2) Needle Free Valve(s) (1) Catheter Securement Device (1) Patient Information Pack (1) Patient Chart Sticker

CATALOG #	NURSING KIT	UOM
MR82013105	3F X 60CM VASCU-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82014105	4F X 60CM VASCU-PICC® VALVED, SINGLE LUMEN	5/BOX
MR82014205	4F X 60CM VASCU-PICC® VALVED, DOUBLE LUMEN	5/BOX
MR82015205	5F X 60CM VASCU-PICC® VALVED, DOUBLE LUMEN	5/BOX

NURSING KIT CONTAINS: (1) PICC (1) Peelable Sheath Introducer (1) Safety Scalpel (1) Guidewire (1) 10cc Syringe (1) Tape Measure (1) Safety Introducer Needle (1|2) Needle Free Valve(s) (1) Catheter Securement Device (1) Latex-Free Tourniquet (1) Patient Information Pack (1) Patient Chart Sticker (1) Twisted Wire Stylet



1499 Delp Dr.

Harleysville, PA 19438

P: 215.256.4201 F: 215.256.1787

www.medcompnet.com

PN2561 Rev. A 9/13