



**F  
R  
E  
E  
F  
L  
O  
W**<sup>™</sup>



**Access. Excellence.**

# FREE-FLOW™

Short & Long Term Hemodialysis Catheter

*As innovation continues  
superior products come to life.*

Medcomp® presents its premier step tip dialysis catheter, the Free Flow™. Through the use of leading technology and unprecedented engineering the industry's first true 360° coaxial flow design has come to life, eliminating the need for side holes and reversing of catheter lines due to arterial insufficiency. The lack of sideholes preserves catheter locking solution concentration and minimizes locking solution leaching. This unique design increases the overall patency and performance of the catheter.



**Raulerson IJ**

## Soft Tip Stylet (Short Term)

Integral soft tip reduces the risk of vessel damage.



**Rotating Suture Wing**  
Rotating suture wing provides external anchoring.




**Crack Resistant Luers**  
The ergonomical design ensures greater resistance to cracking or separation.

**Pre-Loaded Stylet**  
Allows for standard or modified Seldinger introduction. Once removed, renders the venous lumen totally open, unlike tapered tipped catheters.

**True 360° Flow**

- Proven Coaxial Design
- Reduces Arterial Insufficiencies
- No Side Holes Reducing Loss of Locking Fluid

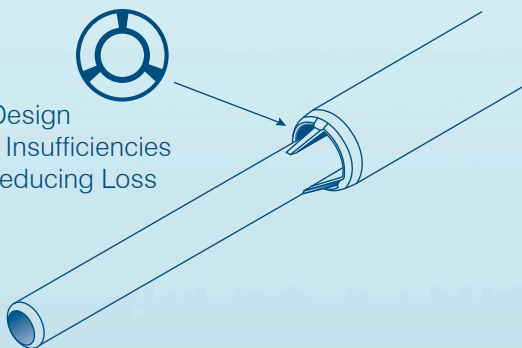

**Open Tips**  
Unobstructed uptake and return reduce pressure and maximize flow.

**FREE FLOW™**

## FREE FLOW™ ORDERING INFORMATION

### LONG TERM HEMODIALYSIS CATHETER

STRAIGHT	Basic Set ( 1 / EA )
15.5FR x 24cm	MR541501
15.5FR x 28cm	MR541502
15.5FR x 32cm	MR541503
15.5FR x 36cm	MR541504
15.5FR x 40cm	MR541505

#### Basic Set

- (1) Catheter w/ Stylet, (1) Grasper Tunneler, (2) Dilators, (1) Scalpel,  
 (1) Valved Tear-Away Sheath, (1) Introducer Needle, (2) Injection Caps  
 (1) Marked .038" Guidewire, (1) Adhesive Dressing,

### SHORT TERM HEMODIALYSIS CATHETER

STRAIGHT	Catheter Only	Basic Set ( 1 / EA )	Full Tray ( 1 / EA )
15.5FR x 15cm	MR541556	MR541551	MR541571
15.5FR x 20cm	MR541557	MR541552	MR541572
15.5FR x 24cm	MR541558	MR541553	MR541573
15.5FR x 28cm	MR541559	MR541554	MR541574

RAULERSON IJ			
15.5FR x 12cm	N/A	MR541561	MR541581
15.5FR x 15cm	N/A	MR541562	MR541582
15.5FR x 20cm	N/A	MR541563	MR541583

#### Basic Set

- (1) Catheter with Stylet, (3) Dilators, (1) Introducer Needle,  
 (1) Marked .038" Guidewire, (1) Adhesive Dressing, (2) Injection Caps

#### Full Tray

- (1) Catheter with Stylet (3) Dilators, (1) Introducer Needle,  
 (1) Marked .038" Guidewire, (1) Adhesive Dressing, (1) Scalpel,  
 (4) Gauze Pads, (1) Chlorhexidine, (3) Safety Hypodermic Needles,  
 (1) Silk Suture, (1) Lidocaine, (1) Hemostat, (2) Syringes,  
 (1) Surgical Gloves, (2) Injection Caps

